pore[®]
enrichment

5 YEAR
WARRANTY

PURERELIEF[®] PRO *EXPRESS*
**FAR INFRARED
VINYL HEATING PAD**

**USER MANUAL
MODEL: PECAIR24-G**

IMPORTANT SAFETY INSTRUCTIONS DO NOT DESTROY

DANGER: TO REDUCE THE RISK OF BURNS, ELECTRIC SHOCK, FIRE, AND PERSONAL INJURY, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS:

1. To ensure the safe and correct use of the product, please read the instructions before use. Please save the instructions properly after reading them.
2. Only use this product on a 110-120V AC circuit. Improper power supply may damage the product.
3. Never pull this pad by the supply cord. Rough pulling may cause damage or electric leakage of the power cord.
4. Do not cut, bend or pinch the cord. Improper use may cause damage, or electric leakage of the power cord.
5. Please check the power cord before each use. Do not use product if the cord is damaged.
6. The heating pad contains electrical circuits. Please do not use any liquid to clean or soak it. When liquid enters the circuit, it may cause a risk of short circuit, product damage, or electric shock.
7. Do not use any conductive metal to prick the heating pad, as there is a circuit inside the pad, which may cause a risk of short circuit or electric shock.
8. Do not use while sleeping.
9. Do not use on an infant or on an animal.
10. Burns can occur regardless of control setting. Check skin under product frequently.

11. Place unit on top of and not under the part of the body that requires heat therapy.
12. Do not use this product with liniment, salve, or ointment preparations as skin burns can result.
13. The product is intended for use in a home environment. The product is not expected to be used together with strong magnetic and high-frequency equipment (such as high-frequency electric knife, microwave oven, short wave equipment, etc.).
14. Do not use clothes to cover this product. Avoid the risk of thermal burns caused by poor heat dissipation conditions.
15. Do not sit on or against the product or crush it. Avoid sharp folds in the fabric.
16. The part of the product that contacts the human body (heating pad) is made of graphene and has passed the tests of ISO 10993-5, ISO 10993-10, and ISO 10993-23 standards. If you have a history of allergy to this product, please avoid direct contact, or stop using it and consult your doctor in case of allergic reaction during contact.
17. This product should not be used in a heat source environment (such as fireplace, indoor heating device). Poor heat dissipation conditions may cause the risk of thermal burns.
18. Do not use in an oxygen-filled atmosphere.
19. Do not use in conjunction with oxygen-rich environments, not intended for use with flammable anesthetics, not intended for use in conjunction with flammable agents.
20. Do not use with generators, power converters, or inverters.
21. Each treatment should not exceed 30 minutes for the same position. Long-term use of the product may cause skin burns.
22. Please keep the product away from children, so as to avoid relevant risks caused by children's play, such as swallowing small parts, damaging the product, or getting electric shock.
23. Avoid children playing with this product. Due to the long power cord, there may be a risk of strangulation or suffocation.
24. This product is not intended for children. Children's skin may cause scalds.

25. Please keep the product away from children. There is a potential risk of suffocation and death when the protective plastic bag is put on the head.
26. During treatment, take good care of children and keep the product away from them.
27. Please keep the product away from pets. Risk of product damage, insulation damage, swallowing small parts, etc. may occur if a pet bites or plays with the product.
28. Please keep the product away from pets during use to avoid device failure and risk of electric shock.
29. Be careful to protect the controller. Do not use the controller carelessly - avoid dropping, hitting, or squeezing it.
30. The controller is not waterproof, please keep it away from liquids.
31. Do not use pins or other metallic objects to fasten this product into place.
32. Do not use the cord as a handle.
33. Unplug when not in use.
34. Loop cord loosely when storing. Tight wrapping can damage cord & internal parts.
35. This product has a polarized plug (one blade is wider than the other). To reduce the risk of electric shock, this plug is intended to fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician to install the proper outlet. Do not attempt to modify the plug or disable this safety feature in any way.
36. The plug is an isolation device that is electrically isolated from the network. Disconnect the electrical connection by unplugging.

Conforms to CAN ICES-003 (B)/NMB-003 (B)

INTRODUCTION

Thank you for purchasing the PureRelief Pro Express Far Infrared Vinyl Heating Pad. We designed this extra-large heating pad to provide temporary soothing heat therapy to large muscle groups in your back, legs, arms, and more.

It features far infrared technology that penetrates muscles more deeply than traditional heating pads with the use of a special graphene fabric cover. Remove the cover, and the heating pad will function as a standard vinyl heating pad.

INDICATION FOR USE

The application of heat therapy is a routinely prescribed therapy in today's medical field. The body's physiological response to heat is dilation of the blood vessels, causing an increase in the blood flow to the area under treatment.

Increased local circulation enhances recovery by flushing away waste products and bringing in fresh blood cells to the treatment area. Heat is exceptionally useful in treating back pain caused by muscle spasms from strain and tension. The pain of arthritic and musculoskeletal conditions can be temporarily alleviated with the use of heat therapy.

The increased blood flow can help relax muscles in spasm and help maintain joint and muscle flexibility. Along with rest, the heating pad will gradually relieve pain by relaxing a muscle in spasm.

Heat treatment is also efficient in relieving pain caused by muscle spasm or inflammation after a day of recreation, gardening, jogging, or household activities.

Intended Use

This heating pad is intended to relieve pain caused by muscle spasm or inflammation after a day of recreation, gardening, jogging, or household activities. It is expected to be used in a homecare environment.

The patient is the intended operator. The use of this product does not require any special skills and training. Anyone who can read this manual can operate it.

Clinical Contraindication

1. DO NOT apply to overly sensitive skin or areas with poor circulation.
2. Individuals with circulation problems should consult with a physician before using this product.
3. This product is not to be used on or by an incapacitated person, a sleeping or unconscious person, a person with diabetes, or a person with poor blood circulation unless carefully attended.
4. Do not use on areas of the skin that are sensitive to touch or heat.

Side Effects

There may be redness or itching on the skin at the treatment site. It can be relieved within 2 hours after applying moisturizer. If discomfort persists, consult your doctor.

Applied Part

This heating pad is classified as a Type BF (body floating) medical device and should only be applied to the skin when used in conjunction with the graphene heating pad cover.

USING THE PURERELIEF PRO EXPRESS HEATING PAD


Preventative Inspection

1. Before each use, check the appearance of the power cord, controller, and heating pad to see if there is obvious mechanical damage, breakage, an abnormal function, and obvious dirt.
2. During electric heating treatment, pay attention to whether the temperature is abnormal and whether there is an abnormal smell.
3. If any abnormality is found during inspection, stop using it immediately and contact your distributor or manufacturer for help.

Instructions

1. Plug the heating pad into a 120V AC electrical outlet.
2. To turn on the PureRelief Pro Express Heating Pad, move the slider button to the Level I heat setting. You can continue adjusting the heat output to Level II and III depending on your desired heat setting. When the heating pad is turned on, the LED light will glow.

I = Low II = Medium III = High

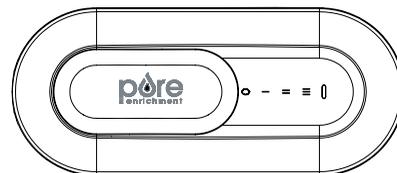3. To turn off, move the slider button to O. The LED light will turn off.

Relationship between LED indication and function:

Function	LED Light Color	Status
Turn On	Green	Waiting Mode
Level I Heat Setting	Green	Low Level
Level II Heating Setting	Green	Medium Level
Level III Heating Setting	Green	High Level

Relationship between heating level and temperature:

Temperature Level	Heat Temperature
Level I Heat Setting	120±20°F (48.88±12°C)
Level II Heating Setting	130±20°F (54.4±12°C)
Level III Heating Setting	140±20°F (60±12°C)

NOTE: Always use the heating pad with the cover on. Do not apply the bare vinyl heating pad directly onto skin.

CLEANING & DISINFECTION MAINTENANCE PROCEDURE

DO NOT DRY CLEAN! DO NOT USE BLEACH! DO NOT USE WRINGER!

NOTE: Long-term use will cause dirt and bacteria on the surface of the graphene cover of the heating pad. Regular cleaning and disinfection is a necessary maintenance procedure to ensure product hygiene.

NOTE: This product is intended for personal use. Please do not share it with others to avoid cross infection.

NOTE: The controller and heating pad of this product contain electrical connection parts, which cannot be cleaned. The cover of the heating pad is a removable and washable part. When the pad cover is dirty, follow the instructions outlined in the "Cleaning and Disinfecting the Cover" section.

Cleaning & Maintaining the Power Cord & Controller

When the surface of the power cord or controller is found dirty by visual inspection, follow the following procedures to clean it:

1. Disconnect the heating pad from the electrical outlet.
2. Use a damp cloth and mild detergent to wipe clean the surface of the power cord and remove stains from the controller.
3. Use a damp cloth to remove any soapy residue.
4. For stubborn stains and to disinfect, soak a clean cloth with 70% Isopropyl Alcohol and wipe clean.

NOTE: The graphene cover material is responsible for the far infrared capability. Over time, machine washing and excess wear and use may reduce the cover's effectiveness.

Cleaning and Disinfecting the Cover

If the surface of the cover cloth of the heating pad is found dirty by visual inspection, or it has been used continuously for more than one month, the following procedures shall be followed for cleaning and disinfection:

1. Disconnect the PureRelief Pro Express Heating Pad from the electrical outlet.

2. Remove the graphene cover from the vinyl heating pad.
3. Machine wash the graphene cover in cold water on the gentle cycle, or hand wash. Line dry or tumble dry on high heat for five minutes. If using an automatic dryer, remove while damp and drape to finish drying.
4. Do not iron the graphene cover.
5. Do not wash the vinyl heating pad or controller.

CALIBRATION, REPAIR & MAINTENANCE

Calibration

This product is a metrological product, which has passed quality inspection before leaving the factory. The expected service life of the product is 5 years. If the product has problems, please replace it with a new one after 5 years.

Repair

NOTE: This product is not authorized for repair or maintenance by any third-party operator. Please contact Pure Enrichment at help@pureenrichment.com if you have issues or concerns using the product.

1. No modification of this equipment is allowed. Improper maintenance or adjustment will reduce the inherent safety and effectiveness of the product.
2. Do not maintain, repair, or adjust the product during use.

Maintenance

The routine maintenance procedures of this product are inspection, cleaning, and disinfection before use. Refer to Pages 8-9 for inspection, cleaning, and disinfection procedures.

COMMON FAULTS AND TROUBLESHOOTING

#	Faults	Reason Analysis of Fault	Solutions	Remarks
1	No response after power on	1. The power socket may be in poor contact 2. Is the power plug plugged in place?	1. Replace a new socket 2. Re-plug	If the problem cannot be solved, please send it to our company for maintenance
2	Push the switch to the corresponding gear, and the indicator does not light up	1. The power plug is not plugged in properly 2. Or plugged in the wrong power supply	1. Re-plug the power supply 2. Identify whether the power supply voltage is AC120V	If the problem cannot be solved, please send it to our company for maintenance
3	No heating	1. It might not be powered on; 2. Equipment hardware failure	1. Identify whether the power indicator is 2. Contact distributor or manufacturer for handling	N/A

PRODUCT SPECIFICATIONS

Power supply	Input: 120V, 60Hz, 50W
Maximum temperature for the heating pad	158°F (70°C) max
Operating Conditions:	41°F (5°C) to 104°F (40 °C), 15-90 % RH, 70-106kPa
Storage and Transportation	-13°F (-25°C) to 158°F (70 °C), 0-90 % RH, 50-106kPa

Recommended maximum use time of single treatment part	30 minutes
Time required for the product to cool down from the highest temperature to 68°F (20°C) between two uses	2 hours
Type of protection against electric shock	Class II
Degree of protection against electric shock	Type BF applied part
Grade of waterproof	IP22
Mode of operation	Continuous operation
Device life	5 years
Note: Not intended to be sterilized. Not for use in an OXYGEN RICH ENVIRONMENT.	
The specification for the fuse	250V T3.15A ; Φ2.7×7.3mm DIP ; E249609/ E221465
The specification for the power cord	C(UL) E309953; SJT 2/C 18AWG(0.824mm2) 221°F (105°C) VW-1 300V
For Heating Pad	
Pad Size	30.5 x 61 cm
Materials and type	Cover: Graphene, 100-180gsm Body: PVC

NORMALIZED SYMBOLS

	Refer to instruction manual/ booklet
	Type BF Applied part
	Class II equipment
	Disposal in accordance with local regulations
IP22	The first number 2: Protected against solid foreign objects of 12.5mm diameter and greater. The second number 2: Protected against vertically falling water drops when enclosure tilted up to 15°.
	Batch Code
	Date of manufacture
	Unique device identifier
	Medical Device

STORAGE AND WASTE DISPOSAL

Storage

It shall be stored according to the product specifications. The product should be stored in a dry and ventilated place that is kept away from children and pets.

NOTE: DO NOT expose the heating pad controller to temperature extremes (below -4°F (-20°C) or over 158°F (70°C) or excessive humidity (>90% RH) as damage could occur to the digital control.

Waste Disposal

After the life of the product expires, the electronic waste of the product shall be disposed according to local regulations. Electronic waste shall not be disposed as domestic waste.

EMC DECLARATION

- The device complies with current specification with regard to electromagnetic compatibility and is suitable for use in all premises, including those designated for private residential purposes. The radio frequency emissions of the device are extremely low and in all probability do not cause any interference with other devices in the proximity.
- It is recommended that you do not place the device on top of or close to other electrical devices, move the unit or connect it to a different socket.
- Radio equipment may affect the operation of this device.

NOTE: The heating pad needs special precautions regarding EMC and needs to be installed and put into service according to the EMC information provided in the ACCOMPANYING DOCUMENTS.

NOTE: Wireless communications equipment such as wireless home network devices, mobile phones, cordless telephones and their base stations, walkie-talkies can effect this equipment and should be kept at least a distance $d=3.3m$ away from the equipment. (Note: As indicated in Table 6 of IEC 60601-1-2:2007 of ME EQUIPMENT, a typical cell phone with a maximum output power of 2W yields $d=3.3m$ at an IMMUNITY LEVEL of 3V/m).

Guidance and manufacturer's declaration - electromagnetic emissions		
The device is intended for use in the electromagnetic environment specified below. The customer or the user assure that it is used in such an environment.		
Emissions Test	Compliance	Electromagnetic environment - guidance
RF emissions CISPR11	Group 1	The device use RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR11	Class B	The device is suitable for use in all establishments other than domestic and those directly connected to the public low voltage power supply network that supplies buildings used for domestic purposes.
Harmonic emissions IEC61000-3-2	N/A	
Voltage fluctuations/ Flicker emissions IEC61000-3-3	N/A	

Guidance and manufacturer's declaration — electromagnetic immunity			
The device is intended for use in the electromagnetic environment specified below. The customer or the user of the device should assure that it is used in such an environment.			
Immunity test	IEC 60601 Test level	Compliance level	Electromagnetic environment-guidance
Electrostatic discharge (ESD) IEC61000-4-2	±8kV direct & indirect contact; ±15kV air discharge	±8kV direct & indirect contact; ±15kV air discharge	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%

Electrical fast transient/bur st IEC 61000-4-4	±2 kV for power supply lines	N/A	N/A (For INTERNALLY POWERED ME EQUIPMENT)
Surge IEC 61000-4-5	± 1 kV line(s) to line(s)	N/A	N/A (For INTERNALLY POWERED ME EQUIPMENT)
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	<5% UT (>95% dip in UT) for 0.5 cycle	N/A	N/A (For INTERNALLY POWERED ME EQUIPMENT)
	40% UT (60% dip in UT) for 5 cycles		
	70% UT (30% dip in UT) for 25 cycles		
	<5% UT (>95% dip in UT) for 5 sec		
Power frequency (50Hz/60Hz) magnetic field IEC 61000-4-8	10V/m	10V/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
Note: UT is the a.c. mains voltage prior to application of the test level.			

Guidance and manufacturer's declaration — electromagnetic immunity			
The device is intended for use in the electromagnetic environment specified below. The customer or the user of the device should assure that it is used in such an environment.			
Immunity test	IEC 60601 Test level	Compliance level	Electromagnetic environment-guidance

Radiated RF IEC 61000-4-3	10V/m & table 9	10V/m & table 9	<p>Portable and mobile RF communications equipment should be used no closer to any part of the Device, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.</p> $d = 1.167 \sqrt{P}$ <p>Recommended separation distance 80 MHz to 800 MHz</p> $d = 2.333 \sqrt{P}$ <p>800 MHz to 2.5 GHz</p> <p>Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in metres (m). Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey, should be less than the compliance level in each frequency range.^b Interference may occur in the vicinity of equipment marked with the following symbol:</p> 
<p>NOTE: At 80 MHz and 800 MHz, the higher frequency range applies. NOTE: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people. a</p>			

A. Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the Device is used exceeds the applicable RF compliance level above, the Device should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the Device.

B. Over the frequency range 150 kHz to 80 MHz, field strengths should be less than [V] V/m.

Test specifications for ENCLOSURE PORT IMMUNITY to RF wireless communications equipment (Table 9)						
Test frequency (MHz)	Band a) (MHz)	Service a)	Modulation b)	Maximum power (W)	Distance (m)	Immunity Test Level (V/m)
385	380-390	TETRA 400	Pulse Modulation b) 18Hz	1.8	0.3	27
450	430-470	GMRS 460, FRS 460	FM c) ±5kHz deviation 1kHz sine	2	0.3	28
710 745 780	704-787	LTE Band 13, 17	Pulse modulation b) 217Hz	0.2	0.3	9
810 870 930	800-960	GSM800/900, TETRA 800, iDEN 820, CDMA 850, LTE Band 5	Pulse modulation b) 18Hz	2	0.3	28
1720 1845 1970	1700-1990	GSM1800; CDMA 1900; GSM 1900; DECT; LTE Band 1,3, 4,25; UMTS	Pulse modulation b) 217Hz	2	0.3	28
2450	2400-2570	Bluetooth, WLAN, 802.11 b/g/n, RFID 2450, LTE Band 7	Pulse modulation b) 217Hz	2	0.3	28

5240	5100-5800	WLAN 802.11 a/n	Pulse modulation b) 217Hz	0.2	0.3	9
5500						
5785						

NOTE: If necessary to achieve the IMMUNITY TEST LEVEL, the distance between the transmitting antenna and the ME EQUIPMENT or ME SYSTEM may be reduced to 1 m. The 1 m test distance is permitted by IEC 61000-4-3.

- a) For some services, only the uplink frequencies are included.
- b) The carrier shall be modulated using a 50 % duty cycle square wave signal.
- c) As an alternative to FM modulation, 50 % pulse modulation at 18 Hz may be used because while it does not represent actual modulation, it would be worst case.

WARRANTY

PureRelief Pro Express comes with an industry-leading 5 Year Warranty that begins on the date of purchase.

The warranty applies to the heating pad and necessary parts and labor relating thereto. The warranty does not apply to damage resulting from failure to follow the operating instructions, accidents, abuse, alterations, or disassembly by unauthorized individuals.

ARE YOU 100% SATISFIED?

If you have product questions or concerns don't hesitate to contact us directly at:

Email: help@pureenrichment.com

Phone: +1 (888) 233-0562

(Available Monday-Friday, 8am-5pm PST)

PRODUCT REGISTRATION

To ensure your full warranty and to receive product updates and streamlined customer support, remember to register your product at:

pureenrichment.com/productregistration

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Fabriquant: Bear Down Brands, LLC