



OPERATION MANUAL READ BEFORE USING



# Meet the Entire Brusheez Family




















# WARNINGS AND PRECAUTIONS READ ALL INSTRUCTIONS.

- This product is not a toy. Children should be supervised by an adult while using this product to avoid risk of choking or misuse.
- This product is designed to clean teeth, gums, and the tongue. Use only for its intended use as described in this manual.
- Discontinue use of this product and contact a physician/ dentist if excessive discomfort or pain is experienced.
- 4. This product is not intended for use by persons with reduced physical, sensory, or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction by a person responsible for their safety.
- Discontinue use if this product appears damaged in any way (brush head, handle, etc.). There are no userserviceable parts. Refer to the Warranty section on Page 7 of this manual if your Brusheez no longer works properly or needs repair.
- 6. Children should be supervised to ensure that they do not play with the product.
- Brusheez is a personal care device and is not intended for use by multiple users.
- 8. Keep the protective cover away from small children as it could present a choking hazard.
- Do not use attachments other than those recommended in this manual.
- To avoid damage to the product, do not place the brush head, handle, or brush covers in the dishwasher for cleaning.
- 11. Always switch the appliance off if:
  - It malfunctions while you are using it.
  - Before you start maintenance.
  - You change or remove accessories.
  - You have finished using the appliance.
- 12. SAVE THESE INSTRUCTIONS.

#### INTRODUCTION

Thank you for purchasing Brusheez. We carefully designed this kids' toothbrush set to make developing healthy dental hygiene habits fun and easy for your little ones between the ages of 3 and 7. Below is a diagram of all of the components included with this toothbrush set.

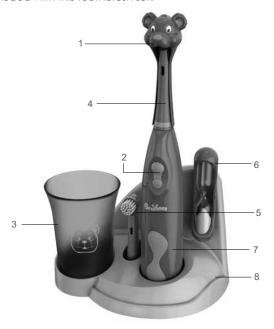

- 1. Protective Cover
- 2. On/Off Switch
- 3. Rinse Cup
- 4. Brush Head

- 5. Extra Brush Head
- 6. Sand Timer
- 7. Brush Handle
- 8. Base

# **USING BRUSHEEZ**

- Fill the rinse cup (3) with water and thoroughly rinse the mouth with it.
- 2. Place a brush head (4) on the brush handle and secure it in place.
- Wet the bristles on the brush head and apply a small amount of toothpaste.
- 4. Place the toothbrush bristles against the teeth at a slight angle towards the gum line.
- 5. Press the On/Off switch to begin brushing. Always insert Brusheez into your mouth before turning it on to prevent the toothpaste from splattering.
- 6. Turn over the sand timer (6) to begin timing the brush cycle. The sand timer is designed to pour the last grain of sand at approximately two minutes. This is the minimum time to brush your teeth thoroughly.
- Move the brush along teeth and gums, applying gentle pressure. Be careful not to press too hard. Let the bristles do the work for you.
- 8. After brushing, switch off the brush and rinse the bristles under the tap. After that, put the protective cover (1) on the brush and place the brush back in the base (8).

# **Weekly Brushing Chart**

A brushing chart can make building healthy oral hygiene habits fun and rewarding. We have included a weekly chart your child can fill in as they learn how to brush their teeth on their own. You can download additional free charts at Brusheez.com.

**NOTE:** For optimal brushing, we advise you to change the brush head every three months or sooner if the brush hairs begin to widen.


The brush head is easy to change by pulling it off the toothbrush. The new brush head can be attached very easily by pushing it on the toothbrush. New brush heads can be purchased at Brusheez.com.

# **USING BRUSHEEZ (CONT'D)**

WARNING: The first few times your child brushes with an electric toothbrush, he or she may experience sensitivity. This will stop within the first week once the gums adjust to the electric toothbrush. If the sensitivity persists longer than a week or returns later after a period of time has passed, please contact your dentist.

#### **CHANGING THE BATTERIES**

Brusheez is powered by two AA batteries (1.5V). For optimal operation, we recommend that you change the batteries once the motor begins to rotate slower than it normally does. Remove the old batteries by opening the battery compartment at the bottom of the brush. Place the new batteries in the battery compartment according to the +/-guides printed inside the unit. Close the lid securely.

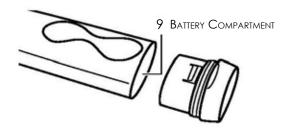

#### **CLEANING & MAINTENANCE**

- Rinse the brush head and bristles after each use. Let the brush head air dry only.
- Remove the brush head once a week and clean the connection between the brush head and the handle.
- Clean the rest of the handle periodically with mild soap and a moist cloth.
- 4. DO NOT use harsh, abrasive cleansers.
- 5. Replace the brush head every 3 months for optimal results.
- 6. When not using Brusheez for an extended period, you should remove the batteries.

#### MOUNTING INSTRUCTIONS

The Brusheez base can be placed on a level surface or mounted on a wall.

# Instructions for wall mounting:

- Before you drill, make sure there are no electrical wires or pipes running along the part of the wall where you plan to install the Brusheez base.
- 2. Drill holes in the wall spaced at the same distance apart as the holes in the back of the base.
- Use the supplied screws and wall plugs to secure the base to the wall.

## Note:

- DO NOT mount the base on the edge of the sink!
- Mounting should be handled by an adult.

#### WARRANTY

Brusheez comes with an industry-leading 5 Year Warranty that begins on the date of purchase.

The warranty applies to the toothbrush set and necessary parts and labor relating thereto. The warranty does not apply to replacement brush heads or damage resulting from failure to follow the operating instructions, accidents, abuse, alterations, or disassembly by unauthorized individuals.

### PRODUCT REGISTRATION

To ensure your full warranty and to receive product updates and streamlined customer support, remember to register your product at:

www.brusheez.com

#### ARE YOU 100% SATISFIED?

If you have product questions or concerns, don't hesitate to contact us directly at:

Email: support@brusheez.com

Phone: (657) 275-3735

(Available Monday-Friday, 8am-5pm PST)



8