

LITTLE
Brushiez[®]
Toddlers' Sonic Toothbrush



Models:

PETODBSH-DNO

PETODBSH-FOX

PETODBSH-KLA

PETODBSH-PDA

PETODBSH-PPY

PETODBSH-UNI

**OPERATION MANUAL
READ BEFORE USING**



Meet the Little Brusheez Family



WARNINGS AND PRECAUTIONS READ ALL INSTRUCTIONS.

1. This product is not a toy. Children should be supervised by an adult while using this product to avoid risk of choking or misuse.
2. This product is designed to clean teeth, gums, and the tongue. Use only for its intended use as described in this manual.
3. Discontinue use of this product and contact a physician/ dentist if excessive discomfort or pain is experienced.
4. This product is not intended for use by persons with reduced physical, sensory, mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction by a person responsible for their safety.
5. Discontinue use if this product appears damaged in any way (brush head, handle, etc.). There are no user-serviceable parts. Refer to the Warranty section on Page 7 of this manual if your Little Brusheez no longer works properly or needs repair.
6. Children should be supervised to ensure that they do not play with the product.
7. Little Brusheez is a personal care device and is not intended for use by multiple users.
8. Keep the protective cover away from small children as it could present a choking hazard.
9. Do not use attachments other than those recommended in this manual.
10. To avoid damage to the product, do not place the brush head, handle, or brush covers in the dishwasher for cleaning.
11. Always switch the appliance off if:
 - It malfunctions while you are using it.
 - Before you start maintenance.
 - You change or remove accessories.
 - You have finished using the appliance.
12. SAVE THESE INSTRUCTIONS.

INTRODUCTION

Thank you for purchasing Little Brusheez. We carefully designed this toddlers' sonic toothbrush set to make developing healthy dental hygiene habits fun and easy for your little ones between the ages of 1 and 3. Below is a diagram of all the components included with this toothbrush set.



- | | |
|---------------------|---------------------|
| 1. Protective Cover | 5. Extra Brush Head |
| 2. Brush Head | 6. Brush Handle |
| 3. Light Timer | 7. Base |
| 4. On/Off Switch | |


USING LITTLE BRUSHEEZ

1. Place a brush head on the brush handle and secure it in place.
2. Wet the bristles on the brush head and apply a small amount of toothpaste.
3. Place the toothbrush bristles against the teeth at a slight angle towards the gum line.
4. Press the On/Off switch to begin brushing. Always insert Little Brusheez into your child's mouth before turning it on to prevent the toothpaste from splattering.
5. The Light Timer will guide your child through the 2-minute brushing process.
6. Move the brush along teeth and gums, applying gentle pressure. Be careful not to press too hard. Let the bristles do the work for you.
7. After brushing, rinse the bristles under the tap. After that, put the protective cover on the brush and place the brush back in the base.

HOW THE LIGHT TIMER WORKS


The 2-Minute Light Timer will guide you and your toddler through the brushing process. It will flash 4 different light colors every 30 seconds, one color per quadrant, to teach toddlers to clean each quadrant. In between each step, the toothbrush will pulse slowly giving your toddler an additional cue that it is time to move on to the next quadrant.

- **Green Light** – Quadrant 1
- **Blue Light** – Quadrant 2
- **Yellow Light** – Quadrant 3
- **Red Light** – Quadrant 4

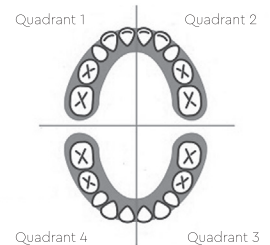

CHANGING THE BRUSH HEAD

The brush head is easy to change by rotating counterclockwise and pulling it off the toothbrush. The new brush head can be attached by placing it on the toothbrush and rotating clockwise

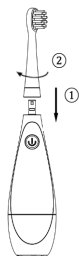


to align in the correct direction. New brush heads can be purchased at Brusheez.com.

NOTE: For optimal brushing, we advise you to change the brush head every 2 months or sooner if the brush hairs begin to widen.

WARNING: The first few times your toddler brushes with a sonic toothbrush, he or she may experience sensitivity. This will stop within the first week once the gums adjust to the sonic toothbrush. If the sensitivity persists longer than a week or returns later after a period of time has passed, please contact your dentist.

CHANGING THE BATTERY

Little Brusheez is powered by a single AAA battery (1.5V). For optimal operation, we recommend that you change the battery once the sonic motor begins to vibrate slower than it normally does. Remove the old battery by opening the battery compartment at the bottom of the brush. Place the new battery in the battery compartment according to the +/- guides printed inside the unit. Close the lid securely.

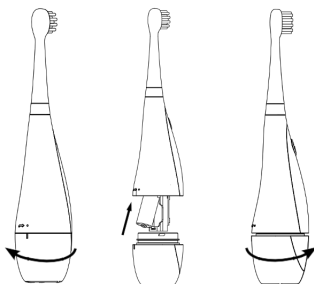

CLEANING & MAINTENANCE

1. Rinse the brush head and bristles after each use. Let the brush head air dry only.
2. Remove the brush head once a week and clean the connection between the brush head and the handle.
3. Clean the rest of the handle periodically with mild soap and a moist cloth.
4. DO NOT use harsh, abrasive cleansers.
5. Replace the brush head every 2 months for optimal results.
6. When not using Little Brusheez for an extended period, you should remove the battery.

WARRANTY

Little Brusheez comes with an industry-leading 5 Year Warranty that begins on the date of purchase.

The warranty applies to the toothbrush set and necessary parts and labor relating thereto. The warranty does not apply to replacement brush heads or damage resulting from failure to follow the operating instructions, accidents, abuse, alterations, or disassembly by unauthorized individuals.

PRODUCT REGISTRATION

To ensure your full warranty and to receive product updates and streamlined customer support, remember to register your product at:

brusheez.com

ARE YOU 100% SATISFIED?

If you have product questions or concerns, don't hesitate to contact us directly at:

Email: support@brusheez.com

Phone: (657) 275-3735

(Available Monday-Friday, 8am-5pm PST)



Santa Ana, CA 92704
brusheez.com

REV: 07.28.22