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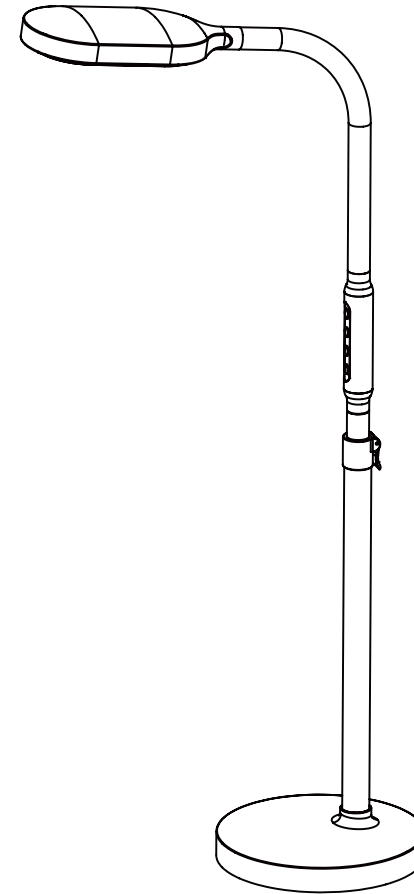
**Model: VT39**

Welcome to the Verilux® family!

We're so excited that you have purchased this HappyLight® Duo 2-in-1 Therapy & Task Floor Lamp. Why? Because light therapy has brightened the lives of millions of people, and you now hold a product in your hands that has the power to brighten your life, too. Pretty amazing stuff!

We want to make sure that you get the most out of your HappyLight® Duo 2-in-1 Therapy & Task Floor Lamp. That's why we stand behind our products and have a three-year limited warranty on this lamp. It's also why we have a super supportive Customer Happiness Team and resources on our website to help you create a light therapy routine that's completely personalized.

Call us at 1-800-786-6850 or visit us at verilux.com if you have any questions as you get started. We're here to help, and we love hearing from customers like you.

Have a bright day!

The Verilux Customer Service Team

## Contents

FCC Notice . . . . .	4
Important Safeguards . . . . .	5-6
Components . . . . .	7
Assembly . . . . .	7
Operation . . . . .	8
Device Pretest . . . . .	8
How to Use . . . . .	9
Establishing a Daily Routine . . . . .	10
Features and Benefits . . . . .	10
Care and Cleaning . . . . .	11
Troubleshooting . . . . .	11
Technical Specifications . . . . .	11
Warranty . . . . .	12

## FCC Notice

Model: HappyLight® Duo 2-in-1 Therapy and Task Floor Lamp (VT39)

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Verilux  
2803 S. Yale St.  
Santa Ana, CA 92704  
verilux.com

The product is protected against unauthorized use. DO NOT modify it without the authorization of the manufacturer, or the user's authority to operate could be voided.

**Note:** This equipment has been tested and found to comply with the limits for a Class B<sup>†</sup> digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

<sup>†</sup> This Class B digital apparatus complies with Canadian ICES-005.

## READ ALL INSTRUCTIONS BEFORE USING

### Important Safeguards

#### **▲ DANGER:**

- To avoid electrocution, do not operate this device near water.

#### **▲ WARNING:**

- Do not use with power supply voltage other than 120 VAC.
- Close supervision is recommended when these lamps are used by or near children, or others who may not understand the need for these precautions.
- To prevent risk of shock or personal injury when cleaning lamp, make sure you have turned off and unplugged the lamp and allowed sufficient time for the lamp to cool.
- Do not shorten or cut the power cord.
- Certain medications (e.g., anti-depressants) and diagnosed mental health conditions may induce light sensitivity, which may result in discomfort, hyperactivity or hypomania. Consult with your healthcare provider before using this product.
- If you have a pre-existing eye condition or are at risk for developing Age Related Macular Degeneration or other eye conditions, consult with your healthcare provider before using this product.
- Some people experience light sensitivity to bright light or upon over-exposure to bright light. Discontinue use if you experience any discomfort.
- Never cover the lamp or place anything on top of it when it is in operation.
- Do not operate this lamp in close proximity to flammable or combustible vapors, such as aerosol spray products, or where oxygen is being administered.

#### **▲ CAUTION:**

- This lamp is for indoor use on a hard surface (e.g., a table or desk) only.
- Do not operate this lamp if it has been damaged in any way.
- Do not disassemble. There are no user-serviceable parts in this lamp.
- This lamp may cause interference with radios, cordless telephones or devices that use a wireless remote control, such as televisions. If interference occurs, move the lamp away from the device, plug the lamp or the device into a different outlet or move the lamp out of the line of sight of the remote control receiver.<sup>†</sup>
- This lamp is not for use with light dimmers, timers, motion detectors or voltage transformers.

<sup>†</sup> This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation."

This Class B digital apparatus complies with Canadian ICES-005.

## Important Safeguards — continued

### **⚠ CAUTION:** — continued

- Use only the power adapter supplied with your lamp. If other power plugs, cords or connectors are used, damage to your lamp or power adapter may occur.
- Unplug lamp during lightning storms or when unused for long periods of time.
- Protect the power cord from being walked on or pinched, particularly at the plug, convenience receptacles and where the power cord jack enters the lamp.
- Do not pull the cable of the AC adapter when unplugging it from the outlet to avoid failure or electric shock.
- Avoid placing the lamp in areas that are dusty, humid/moist, lack ventilation, or are subject to constant vibration.
- Avoid placing this lamp in areas that are exposed to direct sunlight or close to heat-radiating products such as heaters.
- Use a quality surge suppressor on all electronic equipment. Voltage variations and spikes can damage electronic components in any system. A quality suppressor can eliminate the vast majority of failures attributed to surges. Suppressors may be purchased at electronics stores.

**SAVE THESE INSTRUCTIONS**


## Components

### What's Included

- LED lamp
- User manual
- Power adapter
- Allen wrench
- Screw


### Assembly Instructions

1. Find the wire on the gooseneck (A) and insert it into the pole (B), then twist the pole (B) clockwise to connect.
2. Connect the two wires together and insert the pole into the base (C). Tighten the screw using the supplied allen wrench underneath the base (D) after you make sure the touch panel is facing the front.
3. Lift the latch (E) to adjust the height of the lamp pole from 0 to 15 inches.
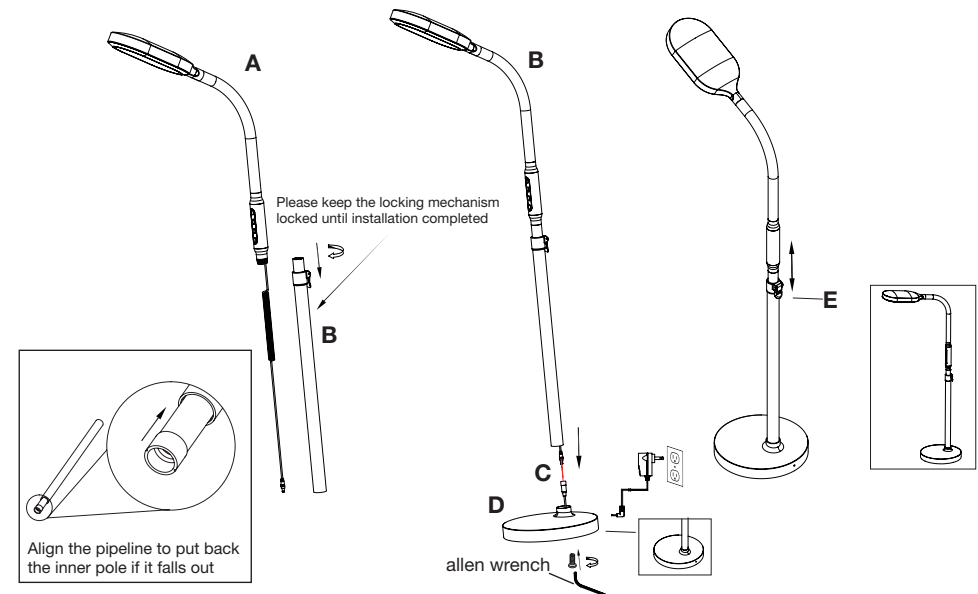
## Operation

### Instructions for Use

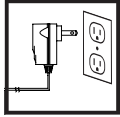**Power Supply:** Plug in the AC adapter into an electrical outlet. Plug the connector of the AC adapter into the HappyLight 2-in-1 Therapy and Task Floor Lamp. (Only use with the supplied AC adapter to avoid damage and fire.)

**On/Off:** To turn on the light, gently touch the on/off touch sensitive control button. (When you turn off the light using the on/off button, it will return to the last setting of brightness and temperature when you turn it on again.)

**Mode:** There are three color temperatures. The initial color temperature is 5000K. To change the temperature simply touch the mode button to change from 5000K (daylight) to 4000K (natural) and then to 3000K (warm).

**Up/Down:** There are five illuminance levels of light intensity on the lamp at each color temperature. Use the Up/Down touch buttons to adjust the illuminance accordingly.

## Device Pretest

To help determine your degree of sensitivity to light, we recommend that you perform this pretest upon the first use of your HappyLight® Duo 2-in-1 Therapy & Task Floor Lamp.

**Step 1** - Turn on the lamp and position it about 18 inches from your face. (Do not look directly into the lamp.)

**Step 2** - Operate the lamp for approximately 10 minutes and turn it off. If the light has caused any discomfort, stop using the lamp and consult a doctor. If you have no discomfort, you may consider gradually increasing the brightness of the light, decreasing your distance from the lamp and increasing your usage time. When you first start using the lamp, you should monitor the degree of any discomfort and accordingly adjust the brightness, distance and time.

**⚠ WARNING:** Certain medications (e.g., anti-depressants) and diagnosed mental health conditions may induce light sensitivity, which may result in discomfort, hyperactivity or hypomania. Consult with your healthcare provider before using this product.

**⚠ WARNING:** If you have a pre-existing eye condition or are at risk for developing Age Related Macular Degeneration or other eye conditions, consult with your healthcare provider before using this product.

## How to Use

### 1 How to Use – Positioning

Sit close to your HappyLight® Duo 2-in-1 Therapy and Task Floor Lamp with your eyes open. While staying oriented toward the light, feel free to engage in other activities such as reading, writing, hobbies and crafts, sewing, needlepoint, applying make-up or working at your computer or desk.

The lamp could be positioned as close as 6 inches but within 24 inches of your face. It is most effective when it is positioned slightly off-center. (See diagram.)

*Note: Do not look directly into the lamp.*

The light from the lamp must be directed at your eyes; therefore, your eyes must be open to achieve the full benefit. *Note: Tinted glasses reduce the amount of light reaching your eyes.*

### 2 How to Use – Customization

Generally speaking, the closer the HappyLight® Duo 2-in-1 Therapy and Task Floor Lamp is to you and the higher the brightness, the less time you will need to use the lamp each day to experience results.

If you do not experience desired results or are feeling lethargic, you can:

- Use the high brightness level
- Sit closer to the lamp
- Spend more time with the lamp

If you experience edginess or are over-stimulated, you can:

- Use the low brightness level
- Sit farther from the lamp
- Spend less time with the lamp
- Take a break from using the lamp

Please note that each person's response to light will vary, so it is important to find the settings, distance and time that are most comfortable for you.

## Establishing a Daily Routine

20 minutes  
to  
1 hour/day\*

Researchers at the National Institute of Health demonstrated that white light, in a range of 2,500 to 10,000 lux (a measurement of light intensity), helped to reverse the symptoms of “Winter Blues.”

We all have a different sensitivity to light, so it is important to find the settings, distance and time

that are comfortable for you. We recommend that you use your HappyLight® Therapy Lamp regularly and at about the same time each day. It is not necessary to sit in front of the lamp continuously. You can divide the time into several mini-sessions. Continue to use the lamp daily to feel an improvement in your mood and activity level.

When you are feeling sluggish and lethargic, increase the brightness of the light, sit closer to the lamp and/or spend more time with the lamp. If you feel a sort of edginess, as though you have consumed too much coffee, it may be time to turn off the lamp.

\* Recommended time is a function of brightness and distance. Individual results may vary. Also, refer to the warnings at the bottom of the Device Pretest section on the previous page.

## Features and Benefits

Your new HappyLight® Duo 2-in-1 Therapy & Task Floor Lamp is a remarkable product designed to offer unparalleled customization for maximum results and comfort. It:

- Delivers full spectrum light to your home, office, dorm room or apartment.
- Cues your body’s own energy enhancers to make you more alert and productive.
- Has up to 10,000 lux<sup>†</sup> to lift spirits, sharpen concentration and improve energy levels.
- Is designed with a compact shape and small footprint to make it easy to place in a convenient location.
- Provides five light brightness levels to assist in gradually increasing light intensity.
- Provides three color temperature options to enhance user comfort.
- Offers simple and easy to read touch controls for power, color temperature, and light brightness.
- Uses long-life, energy efficient LEDs to greatly reduce operating costs over the life of the lamp.

<sup>†</sup>10,000 lux (photopic measurement) at 12".

## Care and Cleaning

Your HappyLight® Duo 2-in-1 Therapy and Task Floor Lamp is made from high-quality materials that will last for many years with minimum care. You may want to periodically clean the lamp using a mild non-abrasive cleaner and soft cloth. When cleaning, make sure you have turned off and unplugged the lamp and allowed sufficient time for it to cool.

## Troubleshooting

If you are experiencing an issue with your HappyLight® Duo 2-in-1 Therapy & Task Floor Lamp, please:

- Make sure the power adapter and plug are inserted fully and securely.
- Make sure there is power to the wall outlet or try another outlet – one that is not connected to a dimmer.

Problem	Check	Solution
Light won't come on.	Outlet end is plugged in properly	Plug the outlet plug into a working outlet.
	Input jack of the power plug is plugged in properly	Be sure it is properly seated in the receptacle on the base.
	Wires in pole and base are properly connected	Be sure wires were properly connected during assembly.

## Technical Specifications

### HappyLight® Duo 2-in-1 Therapy and Task Floor Lamp

**Environment:** Operating Temperature: -20 to 40°C (-4 to 104°F)

**Power:** Input Power to the Adapter: 100-240 VAC, 50/60Hz  
Input Power to the Lamp: DC20V, 2.0A

**Color Temperatures:** Warm: 2700K-3200K  
General Ambient: 3700K-4300K  
Reading/Tasking: 4745K-5311K

**CRI:** >95

**Illumination Intensity:** 10000 Lux

**LED Life Span:** 50000 hours

CETL Listed

RoHS Compliant

Proposition 65 Compliant

## Three-Year Limited Warranty

As the leading innovator of light therapy products for over two decades, Verilux, founded in 1956, backs this HappyLight® Duo 2-in-1 Therapy and Task Floor Lamp with a three-year limited warranty.

**Attention: Once opened, please do not return this product to the store where it was purchased for repair or replacement.**

Many questions can be answered by visiting [verilux.com](http://verilux.com), or you may call our Customer Happiness Team at 1-800-786-6850 during normal business hours.

This limited warranty is provided by Bear Down Brands LLC, 2803 S. Yale St., Santa Ana, CA 92704

Verilux warrants this product to be free from defects in material and workmanship for a period of three years from the date of the original retail purchase from Verilux or an authorized Verilux distributor. **Proof of purchase is required for all warranty claims.** During the limited warranty period, Verilux will replace defective parts of this product when available at no charge to the customer, subject to these limitations: This limited warranty does not include any postage, freight, handling, insurance or delivery fees. This warranty does not cover damage, defect or failure caused by or resulting from accident, external destruction, alteration, modification, abuse, misuse or misapplication of this product.

This warranty does not cover damage to the product resulting from return shipping or handling. Shipping insurance is recommended to help protect your product.

**Return Authorization is required for all returns.** To obtain a Return Authorization, please contact the Verilux Customer Happiness Team.

If, during the first three years of ownership, this product fails to operate properly, it should be returned as specified at [verilux.com/warranty](http://verilux.com/warranty) or as instructed by a Verilux Customer Happiness Specialist.

Due to ongoing improvements, actual product may have slight variations from the product described in this manual.

Please visit our website at: [verilux.com](http://verilux.com) or contact our Customer Happiness Team at

**1-800-786-6850**    [hello@verilux.com](mailto:hello@verilux.com)

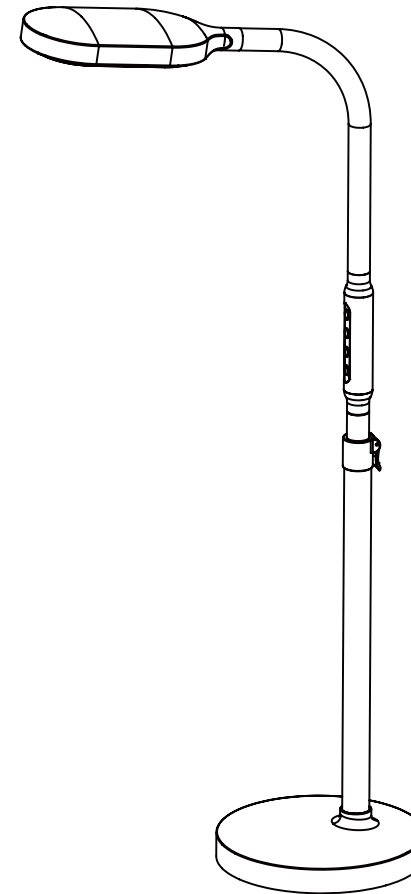
Specialists are available Monday – Friday, 8am – 5pm PST

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