# HappyLight® Alba

LED Light Therapy Lamp  
Lampe de luminothérapie DEL

## Welcome to the Verilux® family!

We're pretty excited that you chose HappyLight® Alba as your new companion! Light therapy has brightened the lives (literally, it's bright!) of millions of people, and you can count on Alba to be there for you when you need it most with healthy, bright light (yeah, we think it's pretty amazing too!).

It's important to us that you get the most out of your experience with HappyLight® Alba. Our website and social media pages are a wealth of information and resources to help you create a light therapy routine that's both effective and personalized. If you find yourself needing a bit more information, don't worry! Our Customer Happiness Team would love to hear from you! And, if anything should go wrong with your light, our three-year limited warranty has got you covered.

Don't hesitate to call us at 1-800-786-6850 or visit our website at verilux.com if you have any questions as you get started. We're here to help you live your best life, with the best light.

Life in a better light,  
Nicholas Harmon  
CEO

#HappyLight

## Contents

<b>Important safeguards</b> .....	<b>4</b>
<b>Setup</b> .....	<b>9</b>
<b>How to Use</b> .....	<b>10</b>
<b>Controls</b> .....	<b>15</b>
<b>Positioning</b> .....	<b>16</b>
<b>Customization</b> .....	<b>17</b>
<b>Device Pretest</b> .....	<b>19</b>
<b>Establishing a Daily Routine</b> .....	<b>20</b>
<b>Features and Benefits</b> .....	<b>22</b>
<b>Care and Cleaning</b> .....	<b>23</b>
<b>Troubleshooting</b> .....	<b>23</b>
<b>Technical Specifications</b> .....	<b>25</b>
<b>Warranty</b> .....	<b>26</b>

## **Important Safeguards**

### **⚠ DANGER :**

---

- To avoid electrocution, do not operate this device near water.

### **⚠ WARNING:**

---

- Close supervision is recommended when these lamps are used by or near children, or others who may not understand the need for these precautions.
- To prevent risk of shock or personal injury when cleaning lamp, make sure you have turned off and unplugged the lamp and allowed sufficient time for the lamp to cool.
- Do not shorten or cut the power cord.

## **Important Safeguards • continued**

- Certain medications (e.g., anti-depressants) and diagnosed mental health conditions may induce light sensitivity, which may result in discomfort, hyperactivity, or hypomania. Consult with your healthcare provider before using this product.
- If you have a pre-existing eye condition or are at risk for developing Age Related Macular Degeneration or other eye conditions, consult with your healthcare provider before using this product.
- Some people experience light sensitivity to bright light or upon over-exposure to bright light. Discontinue use if you experience any discomfort.

### **⚠ CAUTION:**

---

- This lamp is for indoor use on a hard surface (e.g., a table or desk) only.
- Never cover the lamp or place anything on top of it when it is in operation.
- Do not operate this lamp if it has been damaged in any way.

## Important Safeguards • continued

- Do not disassemble. There are no user-serviceable parts in this lamp.
- Use only the power cord adapter supplied with your lamp. If other power plugs, cords, or connections are used, damage to your lamp or power adapter may occur.
- This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Verilux  
14505 Astronautics Lane  
Huntington Beach, CA 92647  
[www.verilux.com](http://www.verilux.com)

- The product is protected against unauthorized use. DO NOT modify it without the authorization of the manufacturer, or the user's authority to operate could be voided.

## Important Safeguards • continued

**Note:** This equipment has been tested and found to comply with the limits for a Class B<sup>†</sup> digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

## Important Safeguards • continued

- This lamp is not for use with light dimmers, timers, motion detectors, or voltage transformers.
- Use a quality surge suppressor on all electronic equipment. Voltage variations and spikes can damage electronic components in any system. A quality suppressor can eliminate the vast majority of failures attributed to surges. Suppressors may be purchased at electronics stores.

**SAVE THESE INSTRUCTIONS**

## Setup and Take Device Pretest

### 1) Setup - Plug In

- Place lamp on a flat, level surface.
- Insert power adapter plug into the adapter socket-as shown on rear of the lamp.
- Plug adapter into a standard 120 VAC electrical outlet.
- Turn on the lamp with the power button located on the top of the lamp. The button will click as you press down.

### 2) Take Device Pretest

Before you begin using your HappyLight® Alba, please take the Device Pretest on page 19 and read through the manual for helpful information and safety precautions.

## How to Use - Controls

Your HappyLight® Alba comes with touch controls, making your light therapy experience easy and fun!

### On/Off

The main power button is located on the top of Alba. The light must be turned on before you can change specific settings. When you turn off your light, your selections are automatically saved for the next time you use your light. How cool is that?

There are three controls that allow you to change your settings and customize your HappyLight® experience:

## How to Use - Controls • continued

### HappyHue™ Color Selector

Use the left control (oriented as if you are looking at the back of the light) to choose between 3 different hues of white light. Each time you touch the control, it will change color.<sup>1</sup>

<sup>1</sup> Color temperature options are warm white (3,500K), neutral white (4,150K), and daylight (5,000K)

## How to Use - Controls • continued

### Brightness

Touch the right control (oriented as if you are looking at the back of the light) to choose between 4 different brightness levels to find the one most comfortable for you.<sup>2</sup> If the brightest setting is uncomfortable, you can start with a lower setting and work your way up to the full 10,000 lux over multiple sessions. Increasing the brightness when in use will allow you to get your daily dose of healthy, bright light more quickly (which means shorter light therapy sessions!). However, if a lower setting is more comfortable for you, that's OK! Just increase your time with the light and you will still get the same benefits!

<sup>2</sup> Brightness level options are 2,500 lux, 5,000 lux, 7,500 lux, and 10,000 lux.

## How to Use - Controls • continued

### Countdown Timer

The countdown timer helps you keep track of the length of your light therapy experience (which makes it pretty convenient). You can set the timer while you do other things (like reading a book, checking email or just enjoying a cup of coffee or tea), and your light will turn off on its own once your light therapy session is complete.

### One Hour Quick Set

To set the timer for a one-hour session, quickly touch the timer control. The illuminated lights will pulse once, indicating that a one-hour session has been set.

## How to Use - Controls • continued

### Customizing the Length of Your Session

You can customize the length of your light therapy session in 10 minute increments.

**Step 1** - Touch and hold your finger on the timer control located on the base of the light for 2 seconds. The illuminated increment lights will pulse twice indicating that programming mode has begun.

**Step 2** - Continue to hold your finger on the control and select the length of your light therapy session. The amount of time is indicated for each light.

**Step 3** - Remove your finger when the desired time has been illuminated. The illuminated lights will pulse three times to recognize the chosen session time and indicate that the programming mode has been completed.

## How to Use - Controls • continued

**Note:** If you touch the control again before the completion of the three pulses, you will be able to continue to modify the session time.

### Stopping the Timer

Touch and hold the timer control until the illuminated lights turn off. If you would also like to turn off the light, press the on/off button at the top of the lamp to turn the power off completely.

## How to Use - Positioning

Sit close to your HappyLight® Alba with your eyes open. While staying oriented toward the light, feel free to engage in other activities that you would normally do. These can include reading, working at your computer or desk, using your phone or tablet, applying make-up, or even hobbies and crafts! Really, there is no limit to what you can do - get your dose of bright light for the day and do your normal routine (multi-tasking #winning)!

You can place HappyLight® Alba between 6 and 24 inches of your face. The time of use and distance are correlated - so the closer you are, the less time you need in your light therapy session and vice versa when it is placed further away (so, more time when it is further from you!). It is most effective when positioned slightly off-center (check out the diagram).

**Important note:** Do not look directly into the lamp.

## How to Use - Customization

Since the light does need to be oriented towards your face (but not staring into it directly), your eyes need to be open to get the full benefit. **Note:** Tinted glasses reduce the amount of light reaching your eyes and would result in the need for a longer session.

We are all different and each person's response to light will vary. It's important to find the distance, time, and color that are most comfortable for you. How close you place HappyLight® Alba to you affects the length of time using the light. The closer you are, the higher the brightness and the less time you will need to use HappyLight® Alba each day to experience the benefits!

If you do not experience desired results, or are feeling lethargic, here are a few things you can try:

- Try using a higher brightness level
- Sit closer to the lamp
- Increase the amount of time you are using the lamp

## How to Use - Customization • continued

If you experience edginess or feel over-stimulated, try some of these options:

- Use a lower brightness level
- Sit farther from the lamp
- Decrease the amount of time you are using the lamp
- Take a break from using the lamp

Don't forget that you have a few amazing options to customize your experience! Alba's HappyHue™ Color Selector lets you adjust the temperature of your light, and the countdown timer helps you easily set your desired session time!

Try keeping a journal of how you feel when using HappyLight® Alba to find the settings that work best for you!

## How to Use - Device Pretest

To help determine your degree of sensitivity to light, we recommend that you perform this pretest upon the first use of your HappyLight® Alba Therapy Lamp.

**Step 1** - Turn on the lamp and position it about 18 inches from your face (remember, don't look directly into the lamp).

**Step 2** - Operate the lamp for approximately 10 minutes and turn it off (you can use the countdown timer to set it for 10 minutes to make things easier!). If the light has caused any discomfort, stop using the lamp and consult a doctor. If you have no discomfort, you may consider gradually increasing the brightness of the light, decreasing the distance from the lamp, and/or increasing your usage time. When you first start using the lamp, you should monitor the degree of any discomfort and accordingly adjust the brightness, HappyHue, distance, and time. Keeping a journal or log of how you feel when using HappyLight® Alba during the first few weeks will help you find the settings, distance, and time that work best for you!

**⚠ WARNING:** Certain medications (e.g. antidepressants) and diagnosed mental health conditions may induce light sensitivity, which may result in discomfort, hyperactivity, or hypomania. Consult with your healthcare provider before using this product.


## How to Use - Device Pretest • continued

**▲ WARNING:** If you have a pre-existing eye condition or are at risk for developing Age Related Macular Degeneration or other eye conditions, consult with your healthcare provider before using this product.

## Establishing a Daily Routine

Establishing your daily routine with HappyLight® Alba will help you to maximize the benefits of light therapy. Clinical research on bright light therapy consistently demonstrates that 2,500 - 10,000 lux (a measurement of the intensity or “brightness” of light) between 20 minutes and a few hours per day help with mood and sleep issues and can reverse the symptoms of “winter blues”. Consistent use at the same time each day has been shown to have the biggest impact and timing earlier in the day has also been shown to be more helpful for sleep cycles.

## Establishing a Daily Routine • continued

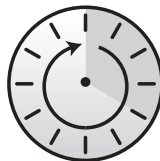20 minutes  
to  
1 hour/day\*

Every person has a different sensitivity to light, so it is important to find the settings, distance, and time that are most comfortable for you. After you have completed the Device Pretest, use your HappyLight® regularly around the same time each day. You could choose to do this all in one session, or divide

your time into several mini-sessions. Use your HappyLight® Alba daily to improve in sleep, boost your mood, increase energy, and enhance your focus.

As you are establishing your routine, if you find yourself feeling sluggish and lethargic, increase the brightness of the lamp, sit closer, and/or spend more time in your session. If you feel a sort of edginess, like you drank too much coffee, it may be time to turn off the lamp.

\*Recommended time is a function of brightness and distance. Individual results may vary. Please also refer to the warnings at the bottom of the Device Pretest section on page 19.

## Features and Benefits

Your new HappyLight® Alba is designed to give you a light therapy experience that's customizable to fit your unique needs. Here are all of the amazing features and benefits you can expect from your new light therapy companion:

- Full spectrum light to brighten up your home, office, or dorm room and cue your body's own energy enhancers to make you more alert, focused, and productive
- 10,000 lux<sup>†</sup> of bright light to improve sleep, boost your mood, increase energy, and enhance your focus
- Adjustable brightness with 4 levels that allow you to choose the intensity that you are most comfortable with
- HappyHue allows you to adjust the temperature of Alba's light from warm white to daylight - find what you like best or change it up depending on how you feel that day
- Thanks to Flicker Elimination Technology® and Optix® Glare Control, you can count on light that is comfortable for your eyes
- Alba's compact shape and small footprint make it easy to find the perfect spot

<sup>†</sup>10,000 lux at 6".

## Care and Cleaning

Your HappyLight® Alba is made from high-quality materials designed to last with minimum care. You may want to periodically clean the lamp using a mild, non-abrasive cleaner and soft cloth. When cleaning, make sure you have turned off and unplugged the lamp and allowed sufficient time for it to cool.

## Troubleshooting

If you are experiencing an issue with your HappyLight® Alba, please:

- Make sure the power adapter and plug are inserted fully and securely.
- Make sure there is power to the wall outlet or try another outlet - one that is not connected to a dimmer.

## Troubleshooting • continued

Problem	Check	Solution
Light won't come on	Have you turned the light on?	The light does not turn on automatically when plugged in. Make sure that you press the power button - it will click (See "How to Use" on page 10)
	Has the plug come loose?	Insert the plug securely into the back of the lamp.
Light Flickers	Has the plug come loose?	Insert plug securely into wall outlet. Insert plug securely into the back of the lamp.
	Is the lamp connected to a dimmer switch?	Connect lamp plug only to wall outlet without dimmer switch

## Technical Specifications

### HappyLight® Alba LED Light Therapy Lamp

**Environment:** Operating Temperature: -10 to 35°C (14 to 95°F)  
Storage Temperature: -20 to 40°C (-4 to 104°F)  
Relative Humidity: 10 to 90% non-condensing

**Power:** Input Power to the Adapter: 100-240 VAC, 0.4A, 50/60Hz  
Input Power to the Lamp: DC20V, 0.8A

**220V Compatible:** For 220V usage, simply attach a plug adapter.  
No converter is needed.

**ETL Listed in Canada and the USA**

## Three-Year Limited Warranty

As the leading innovator of light therapy products for over two decades, Verilux, founded in 1956, backs this HappyLight® Therapy Lamp with a three-year limited warranty.

**Attention: Once opened, please do not return this product to the store where it was purchased for repair or replacement.**

Many questions can be answered by visiting [verilux.com](http://verilux.com), or you may call our Customer Happiness Team at 1-800-786-6850 during normal business hours.

This limited warranty is provided by Verilux, Inc., Huntington Beach, CA 92647.

Verilux warrants this product to be free from defects in material and workmanship for a period of three years from the date of the original retail purchase from Verilux or an authorized Verilux distributor. **Proof of purchase is required for all warranty claims.** During the limited warranty period, Verilux, Inc. will replace defective parts of this product when available at no charge to the customer, subject to these limitations: This limited warranty does not include any postage, freight, handling, insurance, or delivery fees. This warranty does not cover damage, defect or failure caused by, or resulting from, accident, external destruction, alteration, modification, abuse, misuse, or misapplication of this product.

## Three-Year Limited Warranty • continued

This warranty does not cover damage to the product resulting from return shipping or handling. Shipping insurance is recommended to help protect your product.

**Return Authorization is required for all returns.** To obtain a Return Authorization, please contact the Verilux Customer Happiness Team.

If, during the first three years of ownership, this product fails to operate properly, it should be returned as specified at [verilux.com/warranty](http://verilux.com/warranty) or as instructed by a Verilux Customer Happiness Specialist.

Due to ongoing improvements, actual product may have slight variations from the product described in this manual.

## **Three-Year Limited Warranty • continued**

Please visit our website at: **verilux.com** or contact our **Customer Happiness Team** at

**1-800-786-6850**    **hello@verilux.com**

Specialists are available Monday — Friday, 9am — 5pm EST

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La garantie limitée est offerte par Verilux, Inc., Huntington Beach, CA 92647.