
verilux®

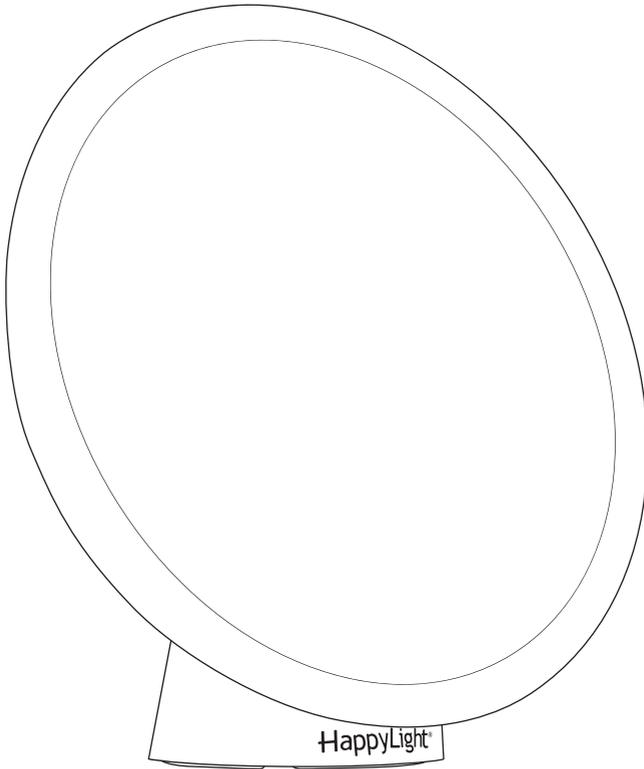HappyLight® halo

Cordless LED Light Therapy Lamp

Model: VT53

Welcome to the Verilux Family!

We're so excited that you have purchased this HappyLight® Therapy Lamp. Why? Because light therapy has brightened the lives of millions of people, and you now hold a product in your hands that has the power to brighten your life, too. Pretty amazing stuff!

We want to make sure that you get the most out of your HappyLight® Therapy Lamp. That's why we stand behind our products and have a three-year limited warranty on this lamp. It's also why we have a super supportive Customer Happiness Team and resources on our website to help you create a light therapy routine that's completely personalized.

Call us at 1-800-786-6850 or visit us at verilux.com if you have any questions as you get started. We're here to help, and we love hearing from customers like you.

Have a bright day!

The Verilux Customer Service Team

Contents

FCC Notice	4
Important Safeguards	5
Setup	6
Device Pretest	6
How to Use	7-11
Charging	7
Plugged In Mode & Unplugged Mode	7
Controls	8-9
Positioning	10
Customization	11
Establishing a Daily Routine	12
Features and Benefits	12-13
Care and Cleaning	13
Technical Specifications	13
Troubleshooting	14
Warranty	14-15

FCC Notice

Model: HappyLight® Halo Cordless LED Light Therapy Lamp (VT53)

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Verilux
2803 S. Yale St.
Santa Ana, CA 92704
verilux.com

The product is protected against unauthorized use. DO NOT modify it without the authorization of the manufacturer, or the user's authority to operate could be voided.

Note: This equipment has been tested and found to comply with the limits for a Class B[†] digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

[†] This Class B digital apparatus complies with Canadian ICES-005.

READ ALL INSTRUCTIONS BEFORE USING

Important Safeguards

▲ DANGER:

- To avoid electrocution, do not operate this device near water.

▲ WARNING:

- Close supervision is recommended when this lamp is used by or near children, or others who may not understand the need for these precautions.
- To prevent risk of shock or personal injury when cleaning the lamp, make sure you have turned off and unplugged the lamp and allowed sufficient time for the bulb to cool.
- Do not shorten or cut the power cord.
- Certain medications (e.g., anti-depressants) may induce light sensitivity which may result in discomfort or a level of hyperactivity. Consult with your health care provider before using this product.
- If you have a pre-existing eye condition or are at risk of developing Age-Related Macular Degeneration or other eye conditions, consult with your healthcare provider before using this product.
- Some people experience light sensitivity to bright light or upon over-exposure to bright light. Discontinue use if you experience any discomfort.

▲ CAUTION:


- This lamp is for indoor use on a hard surface (e.g., a table or desk) only.
- Never cover the lamp or place anything on top of it when it is in operation.
- Do not operate this lamp if it has been damaged in any way.
- Do not disassemble. There are no user-serviceable parts in this lamp.
- Use only the power cord adapter supplied with your lamp. If other power plugs, cords or connections are used, damage to your lamp or power adapter may occur.
- This lamp is not for use with light dimmers, timers, motion detectors, or voltage transformers.
- Use a quality surge suppressor on all electronic equipment. Voltage variations and spikes can damage electronic components in any system. A quality suppressor can eliminate the vast majority of failures attributed to surges. Suppressors may be purchased at electronics stores.

SAVE THESE INSTRUCTIONS

Setup

1. Plug In

- Place lamp on a flat, level surface.
- Insert power adapter plug into the adapter socket – as shown on rear of the lamp.
- Plug adapter into a standard 120V AC electrical outlet.
- Turn on the lamp using the power button located at the back of the lamp as shown.

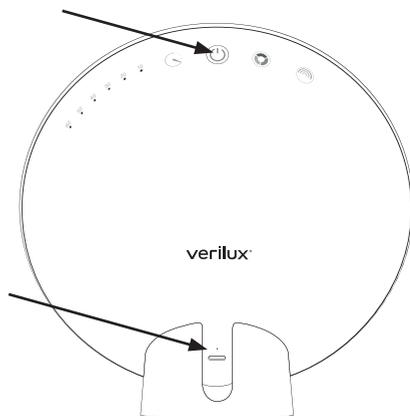Device Pretest

To help determine your degree of sensitivity to light, we recommend that you perform this pretest upon first use of your HappyLight® Halo Therapy Lamp

Step 1 - Turn on the lamp and position it about 18 inches from your face (remember don't look directly into the lamp).

Step 2 - Operate the lamp for approximately 10 minutes and turn it off (you can use the countdown timer to set it for 10 minutes to make things easier!). If the light has caused any discomfort, stop using the lamp and consult a doctor. If you have no discomfort, you may consider gradually increasing the brightness of the light, decreasing the distance from the lamp, and/or increasing your usage time. When you first start using the lamp, you should monitor the degree of any discomfort and accordingly adjust the brightness, HappyHue, distance, and time. Keeping a journal or log of how you feel when using HappyLight® Halo during the first few weeks will help you find the settings, distance, and time that work best for you!

⚠️WARNING: Certain medications (e.g. antidepressants) and diagnosed mental health conditions may induce light sensitivity, which may result in discomfort, hyperactivity, or hypomania. Consult with your healthcare provider before using this product.


⚠️WARNING: If you have a pre-existing eye condition or are at risk of developing Age-Related Macular Degeneration or other eye conditions, consult with your health provider before using this product.

How to Use

Charging

Your HappyLight® Halo will automatically be charging once it is connected to the power adapter and power source. The charging indicator located above the USB-C charging port shows the charging status:


- RED – Charging
 - OFF – Fully Charged
- Charging time from 0% to 100% is approx. 3 hours
 - Cordless usage time varies based on brightness level. The following time estimates are based on a fully charged battery:
 - @25% brightness ~ 2 hours
 - @50% brightness ~ 70 mins
 - @75% brightness ~ 50 mins
 - @100% brightness ~ 40 mins

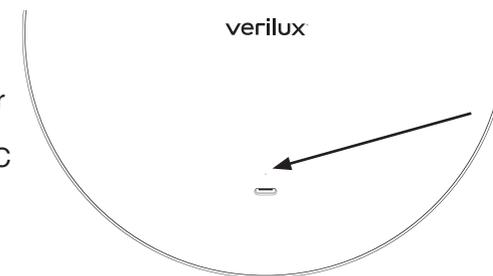How to Use - Plugged-In and Unplugged Mode

Your HappyLight® provides ultimate freedom with a built-in battery that lets you operate the lamp with or without connection to a power source.

- Plugged-In Mode – Simply connect the HappyLight® to the power adapter as shown
- Unplugged Mode – Make sure the HappyLight® is fully charged (Charging Status LED turned OFF), unplug the charging cable and move the HappyLight® to your desired location.

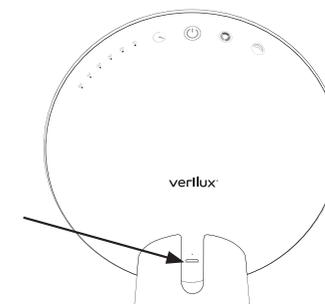How to Use - Controls

Your HappyLight® comes with touch controls, making your light therapy experience easy and fun!

On/Off: The main power button is located on the back of Halo. The light must be turned on before you can change specific settings. When you turn off your light, your selections are automatically saved for the next time you use your light. How cool is that?

There are three controls that allow you to change your settings and customize your HappyLight® experience:

HappyHue™ Color Selector: Choose between 3 different hues of white light. Each time you touch the control, it will change color¹.

¹Color temperature options are warm white (3,500K), neutral white (4,150K), and daylight (5,000K)

Brightness: Choose between 4 different brightness levels to find the one most comfortable for you.² If the brightest setting is uncomfortable, you can start with a lower setting and work your way up to the full 10,000 lux over multiple sessions. Increasing the brightness when in use will allow you to get your daily dose of healthy, bright light more quickly (which means shorter light therapy sessions!). However, if a lower setting is more comfortable for you, that's OK! Just increase your time with the light and you will still get the same benefits!


²Brightness level options are 2,500 lux, 5,000 lux, 7,500 lux, and 10,000 lux.

Countdown Timer: The countdown timer helps you keep track of the length of your light therapy experience (which makes it pretty convenient). You can set the timer while you do other things (like reading a book, checking email or just enjoying a cup of coffee or tea), and your light will turn off on its own once your light therapy session is complete.

One-Hour Quick Set: To set the timer for a one-hour session, quickly touch the timer control. The illuminated lights will pulse once, indicating that a one-hour session has been set.

How to Use - Controls — continued

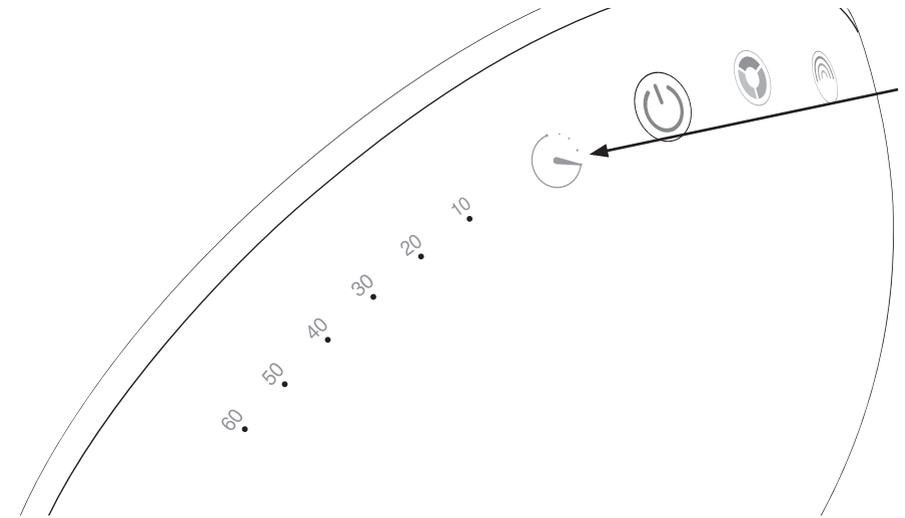Customizing the Length of Your Session

You can customize the length of your light therapy session in 10-minute increments.

Step 1 – Touch and hold your finger on the Countdown Timer button for 2 seconds. The illuminated increment lights will pulse twice indicating that programming mode has begun.

Step 2 – Continue to hold your finger on the control and select the length of your light therapy session. The amount of time is indicated for each light.


Step 3 – Remove your finger when the desired time has been illuminated. The illuminated lights will pulse three times to recognize the chosen session time and indicate that the programming mode has been completed.

Note: If you touch the control again before the completion of the three pulses, you will be able to continue to modify the session time.

Stopping the Timer – Touch and hold the Countdown Timer button until the illuminated lights turn off. If you would also like to turn off the light, press the on/off button on the top of the lamp to turn the power off completely.


How to Use - Positioning

Sit close to your HappyLight® Halo with your eyes open. While staying oriented toward the light, feel free to engage in other activities that you would normally do. These can include reading, working at your computer or desk, using your phone or tablet, applying makeup, or even hobbies and crafts! Really, there is no limit to what you can do – get your dose of bright light for the day and do your normal routine.

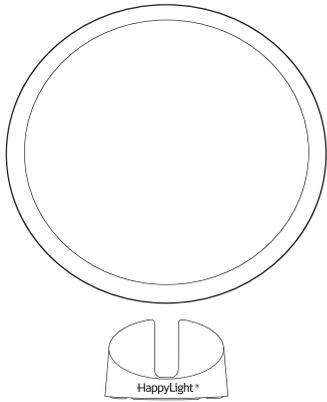Halo's proprietary magnetic and portable base design allows for an adjustable viewing angle simply by rotating it as shown. You can place Halo between 6 and 24 inches from your face. The time of use and distance are correlated - the closer you are, the less time you need in your light therapy session and vice versa when it is placed further away (more time when it is further from you). It is most effective when positioned slightly off-center.

Important Note: Do not look directly into the lamp.

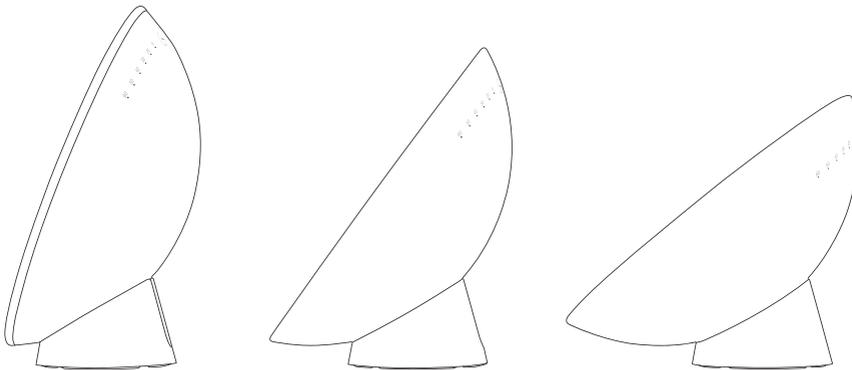How to Use - Customization

Since the light does need to be oriented towards your face (but not staring into it directly), your eyes need to be open to get the full benefit.

Note: Tinted glasses reduce the amount of light reaching your eyes and would result in the need for a longer session.

We are all different and each person's response to light will vary. It's important to find the distance, time, and color that are most comfortable for you. How close you place HappyLight® Halo to you affects the length of time using the light. The closer you are, the higher the brightness and the less time you will need to use HappyLight® Halo each day to experience the benefits!

If you do not experience desired results, or are feeling lethargic, here are a few things you can try:

- Try using a higher brightness level
- Sit closer to the lamp
- Increase the amount of time you are using the lamp

If you experience edginess or feel over-simulated, try some of these options:

- Use a lower brightness level
- Sit farther from the lamp
- Decrease the amount of time you are using the lamp
- Take a break from using the lamp

Don't forget that you have a few amazing options to customize your experience! HappyHue™ Color Selector lets you adjust the temperature of your light, and the countdown timer helps you easily set your desired session time!

Try keeping a journal of how you feel when using HappyLight® Halo to find the settings that work best for you!

Establishing a Daily Routine

20 minutes
to
1 hour/day*

Establishing your daily routine with HappyLight® Halo will help you to maximize the benefits of light therapy. Clinical research on bright light therapy consistently demonstrates that 2,500 – 10,000 lux (a measurement of the intensity

of “brightness of light) between 20 minutes and a few hours per day helps with mood and sleep issues and can reverse the symptoms of “winter blues”. Consistent use at the same time each day has been shown to have the biggest impact and timing earlier in the day has also been shown to be more helpful for sleep cycles.

Every person has a different sensitivity to light, so it is important to find the settings, distance, and time that are most comfortable for you. After you have completed the Device Pretest, use your HappyLight® regularly around the same time each day. You could choose to do this all in one session, or divide your time into several mini-sessions. Use your HappyLight® Halo daily to improve sleep, boost your mood, increase energy, and enhance your focus.

As you are establishing your routine, if you find yourself feeling sluggish and lethargic, increase the brightness of the lamp, sit closer, and/or spend more time in your session. If you feel a sort of edginess, like you drank too much coffee, it may be time to turn off the lamp.

**Recommended time is a function of brightness and distance. Individual results may vary. Please also refer to the warnings at the bottom of the Device Pretest section on page 6.*

Features and Benefits

Your new HappyLight® Halo is designed to give you a light therapy experience that’s customizable to fit your unique needs. Here are all the amazing features and benefits you can expect from your new light therapy companion:

- Full spectrum light to brighten up your home, office, or dorm room and cue your body’s own energy enhancers to make you more alert, focused, and productive
- 10,000 lux[†] of bright light to improve sleep, boost your mood, increase energy, and enhance your focus

Features and Benefits — *continued*

- Adjustable brightness with 4 levels that allow you to choose the intensity that you are most comfortable with
- HappyHue™ allows you to adjust the temperature of Halo’s light from warm white to daylight – find what you like best or change it up depending on how you feel that day
- Thanks to Flicker Elimination Technology and Optix Glare Control, you can count on light that is comfortable for your eyes
- Halo’s compact and small footprint make it easy to find the perfect spot
- Rechargeable battery and cordless operation for portability and convenience
- Magnetic base for flexible light positioning to provide you will the perfect angle of light

[†]10,000 lux at 6”.

Care and Cleaning

Your HappyLight® Halo is made from high-quality materials designed to last with minimum care. You may want to periodically clean the lamp using a mild, non-abrasive cleaner and soft cloth. When cleaning, make sure you have turned off and unplugged the lamp and allowed sufficient time for it to cool.

Technical Specifications

HappyLight® Halo Cordless LED Light Therapy Light

Environment: Operating Temperature: -10 to 35°C (14 to 95°F), Storage Temperature: -20 to 40°C (-4 to 104°F), Relative Humidity: 10 to 90% non-condensing

Power: Input Power to the Adapter: 100-240 V AC, 0.7A max, Input Power to the Lamp: DC12V, 2A max

220V Compatible: For 220V usage, simply attach a plug adapter. No converter is needed.

ETL Listed in Canada and the USA

Troubleshooting

If you are experiencing an issue with your HappyLight® Halo, please:

- Make sure the power adapter and plug are inserted fully and securely.
- Make sure there is power to the wall outlet or try another outlet – one that is not connected to a dimmer.

Problem	Check	Solution
Light won't come on.	Have you turned the light on?	The light does not turn on automatically when plugged in. Make sure that you press the power button – it will click (See “How to Use” on page 7).
	Has the plug come loose?	Insert the plug securely into the back of the lamp.
	Has the light run out of battery?	Make sure the light is connected to the power source and check the charging status LED on the charging level.
Light Flickers	Has the plug come loose?	Insert plug securely into wall outlet. Insert plug securely into the back of the lamp.
	Is the lamp connected to a dimmer switch?	Connect the lamp plug only to wall outlet without dimmer switch.

Three-Year Limited Warranty

As the leading innovator of light therapy products for over two decades, Verilux, founded in 1956, backs this HappyLight® Energy Lamp with a three-year limited warranty.

Attention: Once opened, please do not return this product to the store where it was purchased for repair or replacement.

Many questions can be answered by visiting www.verilux.com, or you may call our Customer Happiness Team at 1-800-786-6850 during normal business hours.

This limited warranty is provided by Bear Down Brands LLC, 2803 S. Yale St., Santa Ana, CA 92704.

Three-Year Limited Warranty — *continued*

Verilux warrants this product to be free from defects in material and workmanship for a period of three years from the date of the original retail purchase from Verilux or an authorized Verilux distributor. **Proof of purchase is required for all warranty claims.** During the limited warranty period, Verilux, will, at its option, repair or replace defective parts of this product at no charge to the customer, subject to these limitations: This limited warranty does not include any postage, freight, handling, insurance, or delivery fees. This warranty does not cover damage, defect, or failure caused by or resulting from accident, external destruction, alteration, modification, abuse, misuse, or misapplication of this product.

This warranty does not cover damage to the product resulting from return shipping or handling. Shipping insurance is recommended to help protect your product.

Return Authorization is required for all returns. To obtain a Return Authorization, please contact the Verilux Customer Happiness Team.

If, during the three years of ownership, this product fails to operate properly, it should be returned as specified at verilux.com/warrantyreplacement or as instructed by a Verilux Customer Happiness Specialist.

Due to ongoing improvements, actual product may have slight variations from the product described in this manual.

**Please visit our website at: verilux.com
or contact our Customer Happiness Team at
1-800-786-6850 hello@verilux.com**

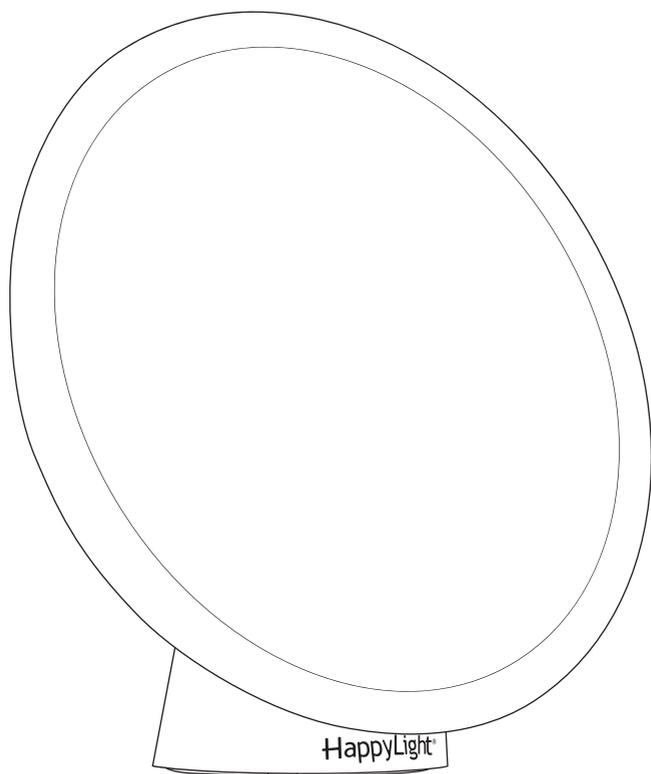
Specialists are available Monday - Friday, 8am-5pm PST

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La garantie limitée est offerte par Bear Down Brands LLC, 2803 S. Yale St., Santa Ana, CA 92704.